




Others



Disposable
Loops / Needles

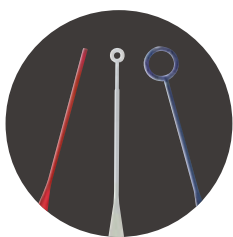

Pipetting Reservoir



Sterilized
Scoop / Spoon


Disposable Loops / Needles

- Disposable inoculation loop and needle made of polystyrene.
- Comfortable and soft-touch handling.
- Available in 3 kinds of package: 1 piece, 10 pieces, 25 pieces.

Note: 10 pieces type and 25 pieces type are packed in a zipper bag.

Real scale

Material: Polystyrene

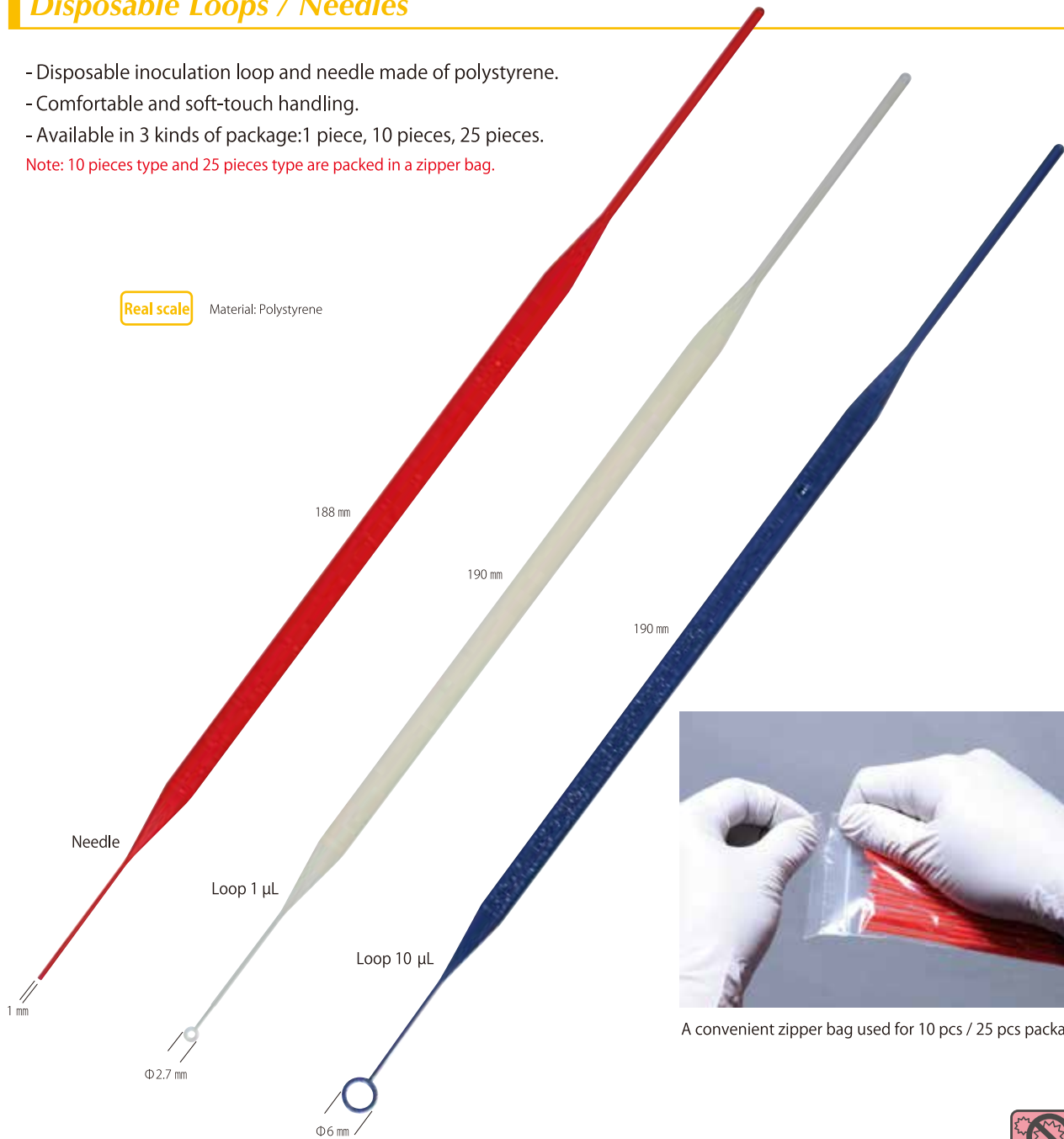

Cat. No.	Item		Unit	Case
171-009-S	Disposable Inoculation Needle, Red	Sterilized	25 pieces × 40 bags	5 units
171-106-S1	Disposable Inoculation Loop, 1 µL, White, Individually wrapped package	Sterilized	250 pieces	5 units
171-109-S10	Disposable Inoculation Loop, 1 µL, White, 10 pcs. package	Sterilized	10 pieces × 100 bags	5 units
171-109-S	Disposable Inoculation Loop, 1 µL, White, 25 pcs. package	Sterilized	25 pieces × 40 bags	5 units
171-406-S1	Disposable Inoculation Loop, 10 µL, Blue, Individually wrapped package	Sterilized	250 pieces	5 units
171-409-S10	Disposable Inoculation Loop, 10 µL, Blue, 10 pcs. package	Sterilized	10 pieces × 100 bags	5 units
171-409-S	Disposable Inoculation Loop, 10 µL, Blue, 25 pcs. package	Sterilized	25 pieces × 40 bags	5 units



Pipetting Reservoir made of polypropylene 20 mL

New Reservoir developed in pursuant of utmost usability.

New feature of "Liquid Sweeper" offers comfortable pipetting and exquisite angle for efficient liquid gathering.



Liquid Sweeper

8 grooves, a very unique feature called "Liquid Sweeper", transfer a liquid drop from the top of the tip quickly.



Drop-on Lid

Special lid for the reservoir to prevent dust or contamination.

(Available only for sterilized type <1791-050-S>)

Note: The lid does not have sealability.

Dimensions: W117mm × D51mm × H6.2mm

■ The most suitable angle for tip touch



The calculated slope angle gathers a liquid almost without its dead volume.

■ Convenient graduation



■ Autoclavable

Note: Minor deformation may occur due to thermal contraction, but it does not affect usage.

■ Easy-to-wash round corners



■ Wide size to cover 12 ch pipette



Cat. No.	Item	Unit	Case
1790-050-30	Pipetting Reservoir, 20 mL, Bulk	30 pieces	10 units
1791-050-S	Pipetting Reservoir, 20 mL, with Lid, Individual package Sterilized	10 pieces	5 units





Disposable Doubled Individually Packaged **Sterilized Scoop / Sterilized Spoon**


- To sample powder or semi-solid material.
- Doubled Individually Packaged is best for using at clean room.
- For easy sample collection available with inside printed scale.
- Flat base design keeps the scoop/spoon upright on the bench.

RNase / DNase / Human-DNA free
 Products with this mark are manufactured in the clean room under our high quality control, and inspected to meet our strict quality standards. RNase/DNase-free certificates can be issued for each lot upon your request.

Pyrogen test
 Results : Endotoxin level less than 0.05 EU/mL
 * according to Internal test

The advantages of doubled individually package
 (e.g., at clean room)

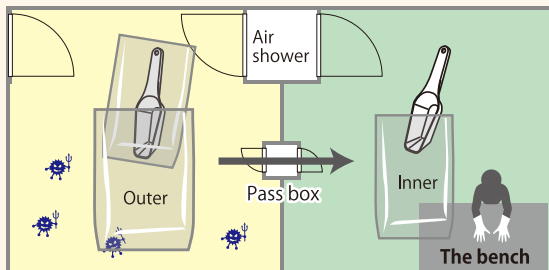

Pre-clean room

Taking outer package out at pre-clean room, then moving to the second room, pass box. Any contaminants or materials on the inner package will be removed in these two rooms.

Clean Room

Taking Scoop/Spoon out from the inner package just before use prevent further contaminants or materials.

Package Material (Outer / Inner): Polyethylene

sterilized SCOOP

100 mL




782-100S

Length Overall: 204.6 mm
 Maximum width of the scoop portion: 59.7 mm
 Height: 31 mm







Material: polypropylene
 (A meet the requirements of FDA)

Measuring scale marked on the inside of the scoop



Self-standing and Flat-bottom design



Cat. No.	Item	Unit
782-100S	100 mL, Sterilized Scoop Sterilized	50 pieces (individual packaging)

sterilized SPOON

1 mL



782-001S

Length Overall: 157 mm
 Maximum width of the spoon portion: 16.9 mm
 Height: 9 mm

5 mL



782-005S

Length Overall: 184.9 mm
 Maximum width of the spoon portion: 24.8 mm
 Height: 9 mm

10 mL



782-010S

Length Overall: 185 mm
 Maximum width of the spoon portion: 29.8 mm
 Height: 18 mm

20 mL



782-020S

Length Overall: 207.7 mm
 Maximum width of the spoon portion: 34.7 mm
 Height: 18 mm

Easy-insert-designed Slim head

Measuring scale marked on the inside of the spoon



Self-standing and Flat-bottom design



Material: polypropylene
 (A meet the requirements of FDA)




Cat. No.	Item	Unit
782-001S	1 mL, Sterilized Spoon Sterilized	100 pieces (individual packaging)
782-005S	5 mL, Sterilized Spoon Sterilized	100 pieces (individual packaging)
782-010S	10 mL, Sterilized Spoon Sterilized	100 pieces (individual packaging)
782-020S	20 mL, Sterilized Spoon Sterilized	100 pieces (individual packaging)

Shape and specifications are subject to change without notice for improvement.

Sterilized Scoop / Sterilized Spoon

